
Peri-implantitis? Management of peri-implantitis without the need for implant removal.

Peri-implantitis? NiTi Brush for peri-implantitis treatment

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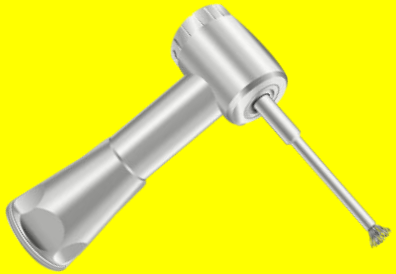
Manage peri-implantitis without implant removal

Treatment without removing the prosthesis

Also for elderly patients & multi-unit bridge cases

NiTi Brush clinical studies (mean age 60–68 yrs, up to 90 yrs)

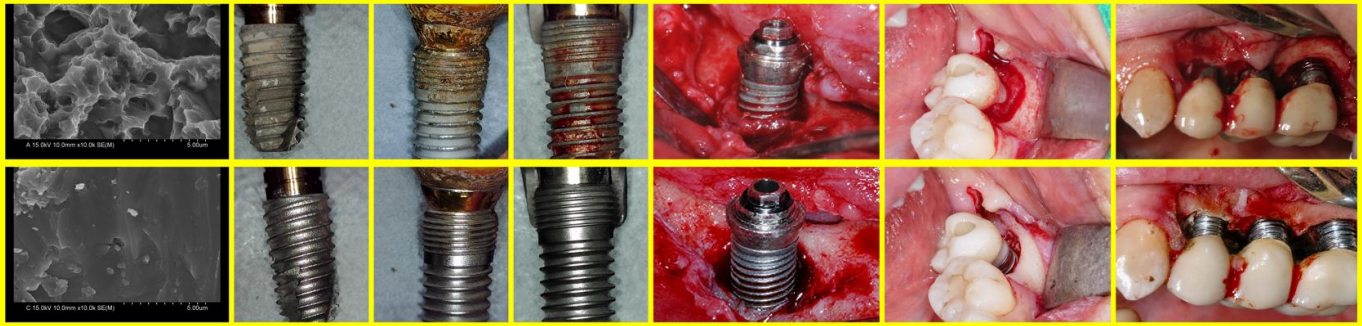HANS Korea

Copyright © 2026 HANS KOREA

www.NiTiBrush.com info@NiTiBrush.com


Peri-implantitis? Management of peri-implantitis without the need for implant removal.

Peri-implantitis? NiTi Brush for peri-implantitis treatment

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Peri-Implantitis? Just brush it away - NiTi Brush

Mount on implant motor RPM 600–1,200 (irrigation required)

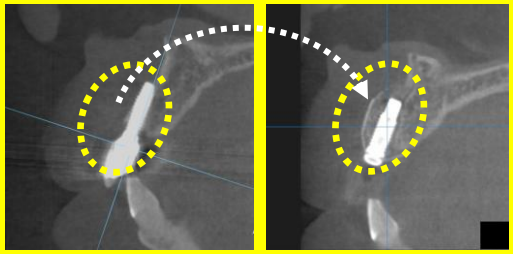Elderly patients & multi-unit bridge cases

NiTi Brush clinical studies (mean age 60–68 yrs, up to 90 yrs)


① Derks J, et al. Reconstructive surgical therapy of peri-implantitis: A multicenter randomized controlled clinical trial. Clin Oral Implants Res 2022;33:921–944.

② Alibegovic L, et al. Bone substitute material in the surgical therapy of peri-implantitis: 3-year outcomes of a randomized controlled trial. Clin Oral Implants Res 2025;36:436–448.

③ Monje A, et al. Surface decontamination on the reconstructive therapy of peri-implantitis: A multicenter randomized clinical trial. Clin Implant Dent Relat Res 2025;27:e70075.

④ Monje A, et al. Extent of implantoplasty in the combined surgical therapy of peri-implantitis: A quasi-randomized clinical trial. Clin Implant Dent Relat Res 2026;28:e70144.

www.nitibrush.com/articles/

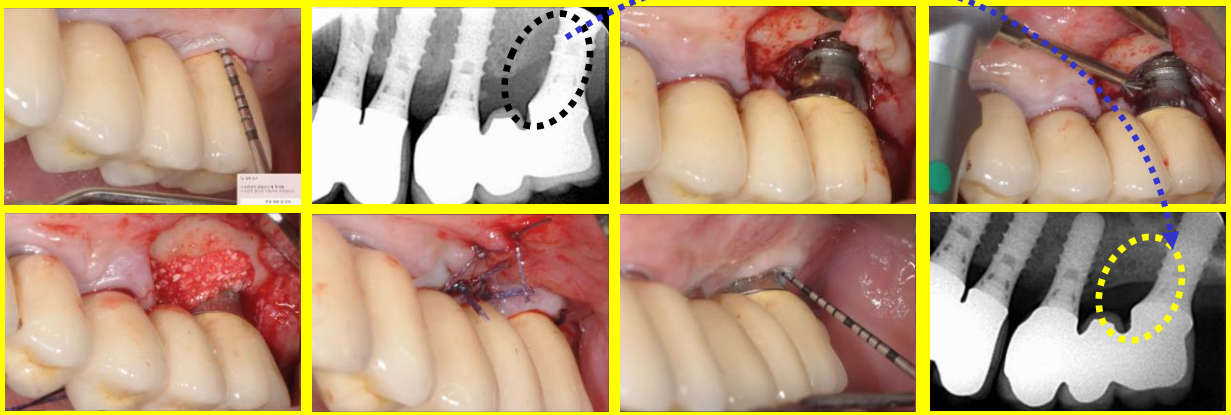HANS Korea

Copyright © 2026 HANS KOREA

www.NiTiBrush.com info@NiTiBrush.com

Peri-implantitis? Management of peri-implantitis without the need for implant removal.

Peri-implantitis? NiTi Brush for peri-implantitis treatment

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
1. Diagnose peri-implantitis

Measure the **depth of pocket**

2. Select a NiTi Brush according to case

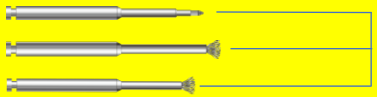Nano for 1 to 2 thread(s) exposure

Pocket for more than 2 threads exposure

Pocket Short for molars

Choice at the clinician's discretion. Pocket: most cases (flap open).
Nano: mild recurrence, without flap opening.

3. Sterilization before use

Before use, sterilize the NiTi Brush for peri-implantitis

Autoclave in sterile pouch; use at once, do not re-store

134 ℃ / 15 minutes or more

Single use · Contraindicated: Ti/Ni allergy

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
4. Rotate NiTi Brush: Surface Cleaning

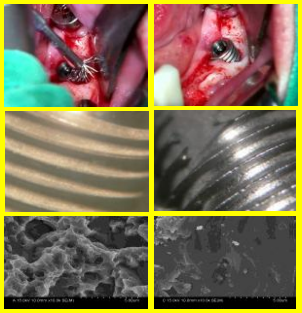Before After

Implant motor

RPM 600 ~1,200

Saline irrigation(external)

Debridement & Decontamination

Using implant motor

Rotating speed is **RPM 600 ~1,200**

Extensive external saline irrigation during brushing

Removes plaque & calculus from the implant surface

Removes necrotic tissue on the implant surface

The cleaned implant surface can be confirmed visually.
Brushing time varies by implant (literature: approx. 2–3 min at 600 rpm).

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
5. Chemical treatment: After brushing (optional)

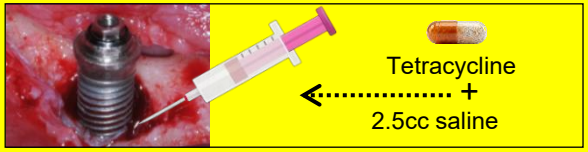Using syringe, apply the solution to the treated site

Wait 3 minutes, then rinse with saline

Open the Tetracycline 250mg capsule

Mix Tetracycline 250mg in 2.5cc saline

Using syringe, apply the solution to the treated site

Wait 3 minutes, then rinse with saline

Literature example. Other agents (e.g., minocycline) have also been reported; saline alone in many studies. Related literature available on request.

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