HANS Korea


www.NiTiBrush.com
info@nitibrush.com

NiTi Brush for peri-implantitis

Make rough surface into machined smooth surface

Remove everything except the original shape of implant



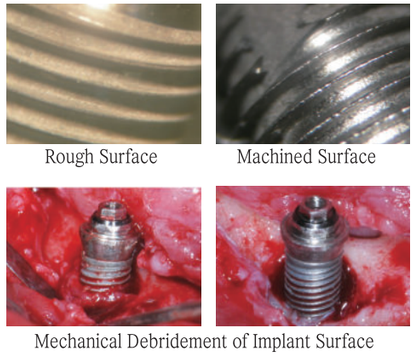


Look for distributors !!!

Rev.002 © 2023 HANS Korea Co., Ltd. All rights reserved.

Selection of NiTi Brush for peri-implantitis

Depend on the case, select one model of NiTi Brush as below.

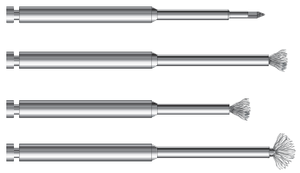

NANO for One thread open

POCKET for 2~3 threads open

POCKET short for molar

OMEGA for deep pocket



Sterile before use(134°C/15min)

Using “implant motor”, **RPM 600~1,200**

Single time use only

Save Implant & Create New Patient

Protocol of the NiTi Brush for peri-implantitis

1. Check up the periimplantitis.




Check up the periimplantitis using probe and possible any solution.

2. Sterile the NiTi Brush.



Before use, sterile the NiTi Brush using sterile pouch & autoclave 134 ℃ / 15 minutes

3. Attach the NiTi Brush to handpiece.

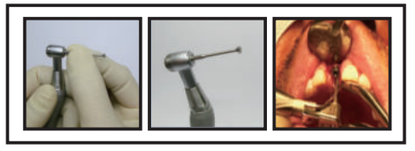


Attach the NiTi Brush to (20:1) implant handpiece
Prepare external saline imigation

4. Brushing – mechanical treatment.



implant motor, RPM 600~1,200
Rotate the NiTiBrush
Remove osseous defect of plague calculus or any other materials
Entensive extenmal sterile saline

5. Check up the mechanical debridement by naked eyes.

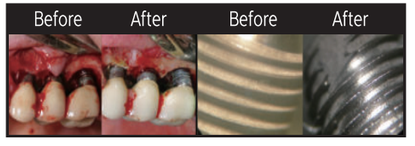

During brushing, you can see
The shade of implant surface changed
into shiny & machined surface